



MedSource
I.V. Infusion

TrueSafe Comfort™
Safety I.V. Catheter



RapidFlash™ Needle

The notched needle, **featured in sizes 20g-24g**, allows for catheter blood flash observation before blood reaches the flash chamber.

Our Safety I. V. Catheter is designed to reduce hit-and-miss insertion attempts by immediately confirming vessel entry, enabling you to stay focused on the insertion site. The TrueSafe Comfort™ is manufactured with PUR material which insures a more comfortable experience for the patient and a longer in-dwelling performance for the catheter.

The design features a notched needle which allows for catheter blood flash observation; incorporated into gauge sizes 20g - 24g. After threading, with the press of a button, the needle instantly retracts into the barrel chamber helping reduce accidental needle sticks and making sharps disposal easier.

Rev. 1

www.gomedsource.com

Item #	Description	UOM
MS-84314	TrueSafe Comfort, 14 G x 1.75"	EA, 50 EA/BX, 4 BX/CS, 200 EA/CS
MS-84316	TrueSafe Comfort, 16G x 1.77"	EA, 50 EA/BX, 4 BX/CS, 200 EA/CS
MS-843182	TrueSafe Comfort, 18G x 1.16"	EA, 50 EA/BX, 4 BX/CS, 200 EA/CS
MS-84318	TrueSafe Comfort, 18G x 1.88"	EA, 50 EA/BX, 4 BX/CS, 200 EA/CS
MS-843202	TrueSafe Comfort, 20G x 1.16"	EA, 50 EA/BX, 4 BX/CS, 200 EA/CS
MS-84320	TrueSafe Comfort, 20G x 1.88"	EA, 50 EA/CS, 4 BX/CS, 200 EA/CS
MS-843203	TrueSafe Comfort, 20G x 1.00"	EA, 50 EA/CS, 4 BX/CS, 200 EA/CS
MS-84322	TrueSafe Comfort, 22G x 1.00"	EA, 50 EA/CS, 4 BX/CS, 200 EA/CS
MS-84324	TrueSafe Comfort, 24G x 0.75"	EA, 50 EA/CS, 4 BX/CS, 200 EA/CS

Features:

- Fully encapsulating needle shield with push button activation
- Notched needle in sizes 20g - 24g
- Contoured housing and ergonomically designed thumb grip
- Gently tapered catheter tip
- Large-volume chamber


Benefits:

- Designed to help health care workers as well as hospital waste handlers from exposure to blood borne pathogens via needle stick injuries.
- Catheter blood observation before blood reaches flash chamber.
- Allows for a variety of user techniques to aid insertion.
- Provides smooth transition from needle to catheter.
- Allows clear visual confirmation through flashback and easy handling.

Specification: TrueSafe™ Safety I.V. Catheter

- Catheter Manufactured With PUR (Poly-Urethane)
- Siliconized Polished Stainless Steel Needle
- RapidFlash™ Notched Needle

- Smooth surface for easy catheter insertion and advancement. More comfortable for patients and a longer in-dwelling performance.
- Sharp needle technology provides both ease of penetration and improved patient comfort during cannulation.
- Available in sizes 20g - 24g.

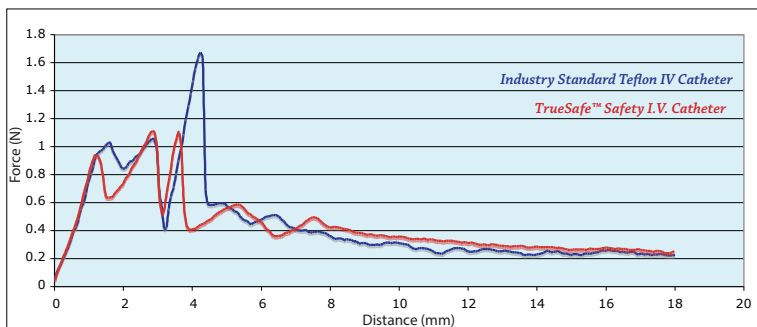

Needle Penetration Comparison

When compared to the industry standard I.V. Catheter, the MedSource Safety I.V. Catheter outperforms in comfort and gives a less painful needle insertion.