



ETCO₂ Sampling Cannulas

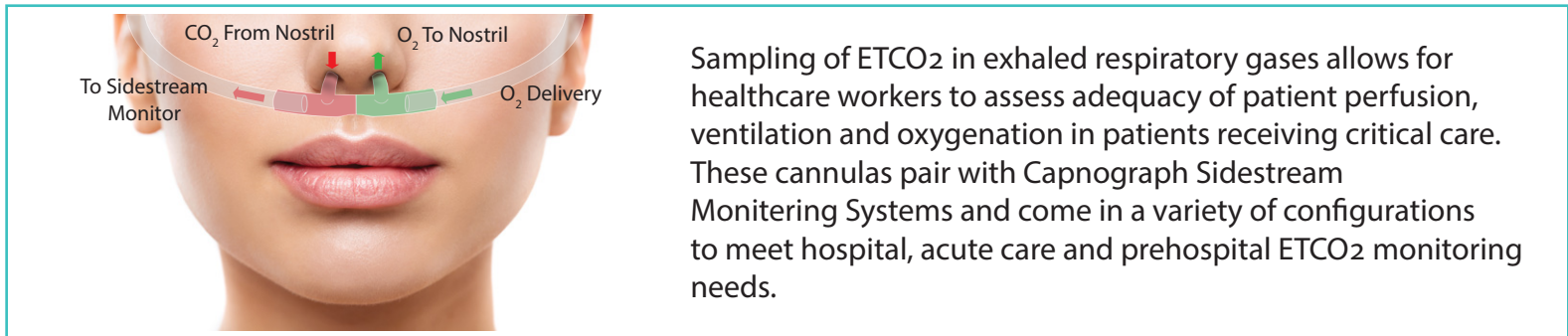

MedSource ETCO₂ Sampling Cannulas allow for a breath-by-breath sampling of ETCO₂ while simultaneously delivering oxygen.

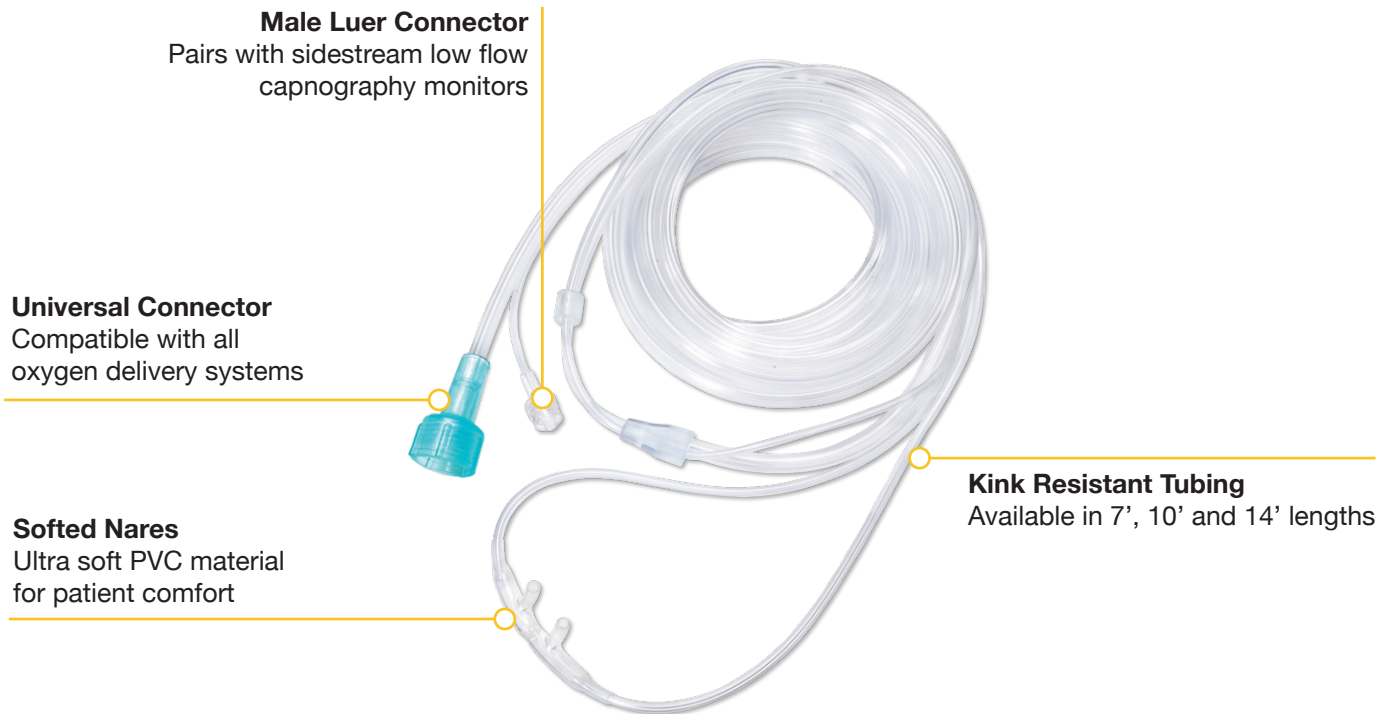

SPECIFICATION

- Dip molded PVC nasal cannula
- Male luer lock
- Not made with natural rubber latex
- PVC material not made with DEHP
- Kink resistant tubing
- Pair with any sidestream monitor (low flow 50ml/min)



www.medsourcelabs.com



Item #	Description	UOM
MS-35400A7	Adult Nasal CO ₂ Sampling Cannula W/O ₂ delivery Luer Conn. 7ft.	25 EA/BX, 200 EA/CS
MS-35400A10	Adult Nasal CO ₂ Sampling Cannula W/O ₂ delivery Luer Conn.10ft.	25 EA/BX, 150 EA/CS
MS-3400A14	Adult Nasal CO ₂ Sampling Cannula W/O ₂ delivery Luer Conn.14ft.	25 EA/BX, 100 EA/CS

ETCO₂ Sampling Cannulas



Dual Nasal/Oral

MedSource ETCO₂ Sampling Cannulas allow for a breath-by-breath sampling of ETCO₂ while simultaneously delivering oxygen.

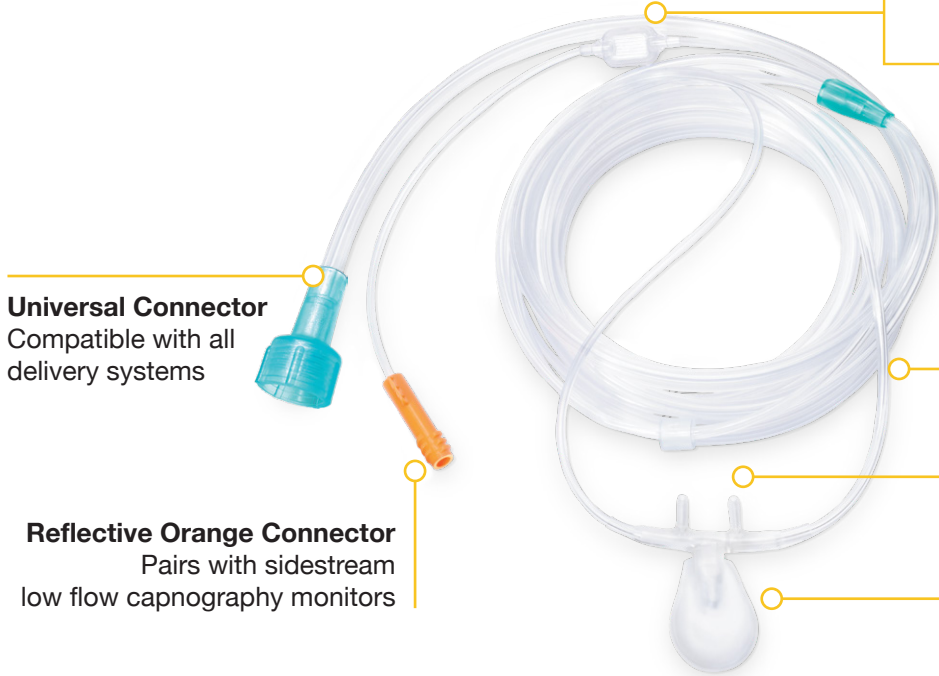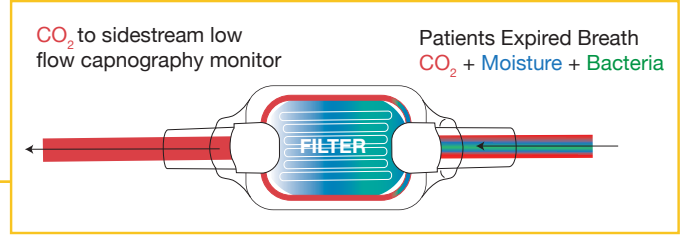

Universal Connector
Compatible with all delivery systems

Reflective Orange Connector
Pairs with sidestream low flow capnography monitors


Button Filter

Combines hydrophobic bacterial and hydrophilic moisture-wicking technology to minimize moisture and bacteria to help ensure an accurate diagnosis

Kink Resistant Tubing

Softed Nares
Ultra soft PVC material for patient comfort

Softed Tear Shaped Scoop
Lifts up and away from patient's lip for extended comfort

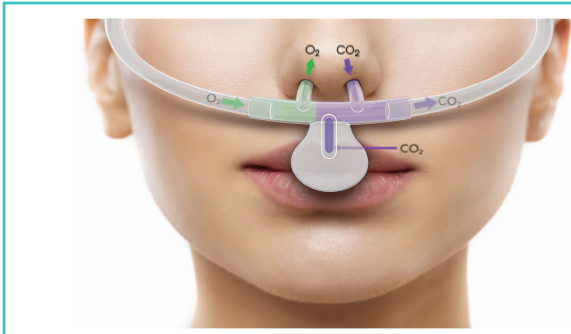

Sampling of ETCO₂ in exhaled respiratory gases allows for healthcare workers to assess adequacy of patient perfusion, ventilation and oxygenation in patients receiving critical care. These cannulas pair with Capnograph Sidestream Monitoring Systems and come in a variety of configurations to meet hospital, acute care and prehospital ETCO₂ monitoring needs.

SPECIFICATIONS

- Dip molded PVC nasal cannula and oral scoop
- Reflective orange connector
- Not made with natural rubber latex
- PVC material not made with DEHP
- 7' of kink resistant tubing
- Pair with sidestream monitor (low flow 50ml/min)



www.medsourcelabs.com

Item #	Description	UOM
MS-35601A	Adult Nasal/Oral Dual CO ₂ Sampling Cannula W/ O ₂ Delivery	200/CS, 25/BX
MS-35601P	Pediatric Nasal/Oral Dual CO ₂ Sampling Cannula W/ O ₂ Delivery	200/CS, 25/BX

ETCO₂ Sampling Cannulas T-Connector



MedSource ETCO₂ Sampling Cannulas allow for a breath-by-breath sampling of ETCO₂.


Button Filter

Minimizes entrance of moisture from sampling ports to monitor

Reflective Orange Connector

Pairs with sidestream low flow capnography monitors

Kink Resistant Tubing

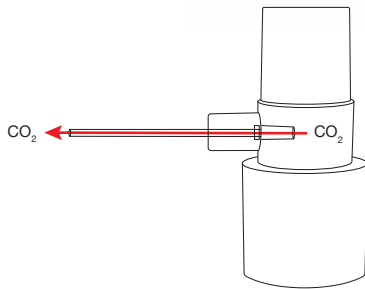

Sampling of ETCO₂ in exhaled respiratory gases allows for healthcare workers to assess adequacy of patient perfusion, ventilation and oxygenation in patients receiving critical care. These cannulas pair with Capnograph Sidestream Monitoring Systems and come in a variety of configurations to meet hospital, acute care and prehospital ETCO₂ monitoring needs.

SPECIFICATION

- Reflective orange connector
- 7' of kink resistant tubing
- Not made with natural rubber latex
- PVC material not made with DEHP
- Pair with sidestream monitor (low flow 50ml/min)



www.medsourcelabs.com

Item #	Description	UOM
MS-35770	Intubated Sampling Line W/T-Connector, Orange Luer W/Reflective Conn., Sampling Scoop, Filter, 7ft.	100/CS, 25/BX