


# NIK Kit

## Needle Decompression Indicator Kit

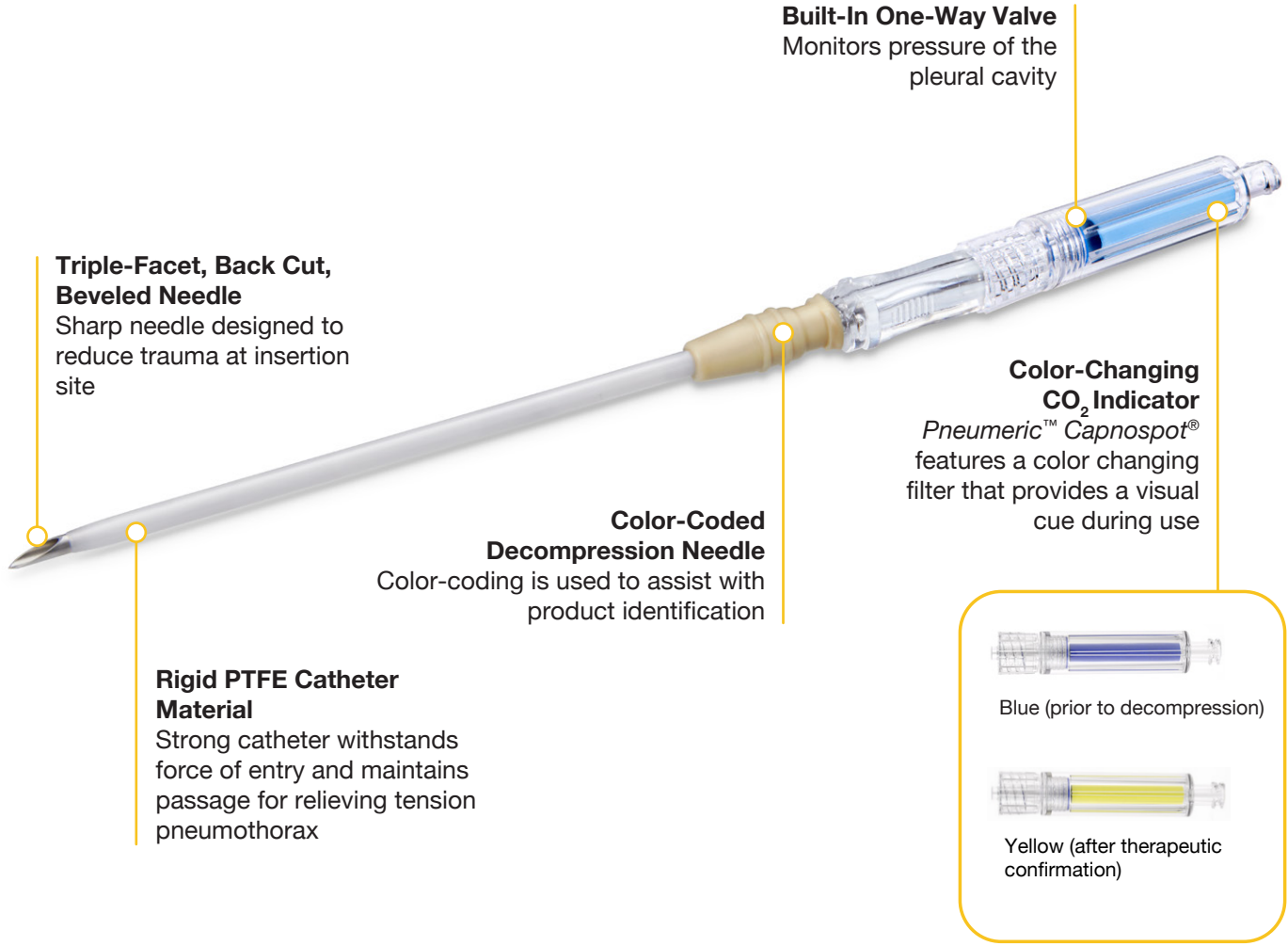

**Triple-Facet, Back Cut, Beveled Needle**

Sharp needle designed to reduce trauma at insertion site

**Built-In One-Way Valve**  
Monitors pressure of the pleural cavity

**Color-Changing CO<sub>2</sub> Indicator**

*Pneumeric™ Capnospot®* features a color changing filter that provides a visual cue during use

**Color-Coded Decompression Needle**

Color-coding is used to assist with product identification

**Rigid PTFE Catheter Material**

Strong catheter withstands force of entry and maintains passage for relieving tension pneumothorax



Blue (prior to decompression)



Yellow (after therapeutic confirmation)



**Buy American Act Compliant**- Manufactured in accordance with federal sourcing requirements for domestic procurement

### SPECIFICATIONS

- Needle Size: 14,12, 10 Gauge x 3.25 inch
- PTFE Catheter material
- EO Sterilization
- Single Use Only
- Not made with natural rubber latex
- Not made with DEHP
- Each kit includes: One Decompression Needle and One *Pneumeric™ Capnospot®*

Item #	Description	UOM
MS-87110	NIK Kit - 10G x 3.25"	EA
MS-87112	NIK Kit - 12G x 3.25"	EA
MS-87114	NIK Kit - 14G x 3.25"	EA
MS-871142	NIK Kit - 14G x 2"	EA
MS-841143	NIK Kit - 14G x 1.5"	EA